

ROLLING STORAGE CART FOR VACUUM CUSHIONS WITH S-HOOKS



Article No. : **32790**

32791 (5pcs)

A. GENERAL PRODUCT INFORMATION

The product referred to in these instructions is a storage system for medical devices used for patient positioning and immobilisation in radiation therapy.

This product may only be used in combination with Orfit immobilisation products. Orfit prohibits the use of unauthorised third-party products in conjunction with its own products.

B. PRODUCT DESCRIPTION

This ROLLING STORAGE CART FOR VACUUM BAGS with S-shaped hooks (5 pieces) is used to store ORFIT vacuum bags.

C. PRECAUTIONS FOR USE

- Always use the correct tools for screwing.
- Make sure that the screws are correctly fastened and no parts are hanging loose
- Wear eye and ear protection and working gloves whenever using power tools.
- Use the breaks on the wheels of the storage cart when not in use.

D. STORAGE

Always store the product in a safe place to prevent it from getting damaged or falling onto other objects.
Store the device between +10°C (50°F) and 40°C (122°F).

E. PROPERTIES

Physical Properties

The following physical properties apply to these units:

Dimensions: D 590 mm x W 1575 mm x H 1600 mm

Weight: 15.5 kg

F. ASSEMBLY INSTRUCTIONS


The storage cart must be assembled by two persons.
Create a safe and clean work area before opening the packaging and handling of all the parts.

G. REQUIRED TOOLS

- Allen keys
- Torqx wrench

H. PARTS LIST

ITEM	QTY	DESCRIPTION
1	2	Carrier plate
2	1	Oval connecting tube
3	2	Carry stand
4	4	Wheel 2x with brake 2x without
5	1	Stainless steel suspension tube
6	2	Cover cap square black
7	6	Cover cap round black
8	4	Broached Socket Head Cap Screw – Metric Taprite M5 x 25
9	2	Broached Socket Head Cap Screw – Metric M8 x 25
10	6	Torx socket countersunk and button head screws – Metric series M10 x 30
11	5	S-hooks

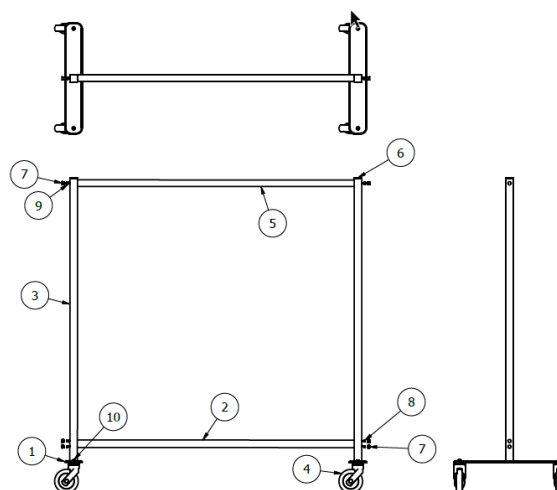




I. ASSEMBLY

1. Attach carrier plate (3) to carry stand (1) using the M10 screws (10). Repeat for the other side.



2. Attach suspension tube (5) to the top of the carry stand (1) using M8 screws (9).

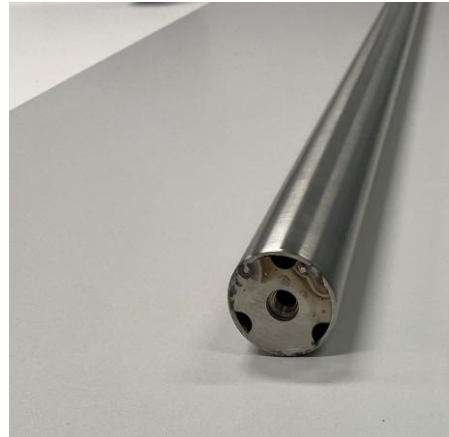

3. Attach connecting tube (2) to the bottom of the carry stand (3) with M5 screws (8). Make sure to properly fix the screws to ensure stability of the cart.



4. Cover openings with cover caps (6,7).





5. Hang S-hooks (11) over the suspension tube (5)

J. MAINTENANCE AND WASTE MANAGEMENT

This product can be cleaned and disinfected by means of soapy water or an isopropanol based disinfectant, applied with a soft cloth. If unsure about the cleaning fluid, do not use. **Never use aerosol sprays, corrosive cleaning agents, solvents or abrasive detergents.**

Periodic checks of the product should be done to insure the parts are not worn and require repair or replacement. **Do not attempt to make repairs yourself.** Contact your distributor if there are any questions or concerns.

The product can be disposed of with household waste.

K. ADDITIONAL INFORMATION

For additional information such as distributor contact information, product brochures, Safety Data Sheets and regulatory information, please visit our website www.orfit.com.

Note:

It is prohibited to make alterations to this text without prior approval from Orfit Industries.



ORFIT INDUSTRIES
Vosveld 9A | B-2110 Wijnegem | Belgium
T (+32) (0)3 326 20 26
welcome@orfit.com

www.orfit.com



Ref. No. 50256
VERSION 2
LAST UPDATE: 08/09/2026
REVISION DATE: 08/09/2028