

THERMOPLASTIC EDGING STRIP

GENERAL PRODUCT INFORMATION

THERMOPLASTIC EDGING STRIP is a thin low-temperature thermoplastic lining material.

It is intended for use as an edging and finishing material in custom-made fabrication applications to provide a smooth finish and improve edge protection.

PRODUCT RANGRE

THERMOPLASTIC EDGING STRIP is available on rolls of 20 mm (width) x 3000 mm (length) x 0.8 mm (thickness).


Art. no.	Colour
360201N.ST1	White

PRECAUTIONS BEFORE USE

- Ensure the workplace is well ventilated.
- Use appropriate tools and ensure safe handling during preparation and assembly.
- Ensure a comfortable and stable working position during use.

WORKING TECHNIQUE


1. Pull the loose end of the edging strip out of the package incision and cut the required length of material using suitable scissors.

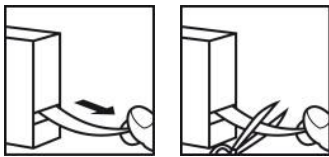


2. Activate the strip in a hot water bath, using a dry heater, or in a hot air oven at a minimum temperature of 65°C (149°F).



3. Wait at least 1 minute until the material becomes transparent. This indicates that the material is ready for moulding.
4. Remove the material from the water bath and briefly dry it with a towel to remove excess water drops.
5. Apply the EDGING STRIP to the edges of the fabricated item after trimming.

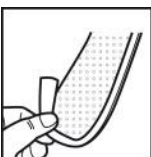

6. Ensure that half of the strip covers the inner surface of the edge and the other half covers the outer surface of the edge.

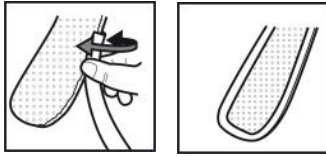

7. EDGING STRIP can be removed after application by carefully reheating the applied edges.

MAINTENANCE

- EDGING STRIP can be cleaned with lukewarm water and liquid soap. Rinse thoroughly and dry completely.

⚠ If you are not sure about the cleaning agent, do not use it.

⚠ Do not use solvents or abrasive detergents.

- EDGING STRIP is not suitable for sterilization in an autoclave.

WASTE MANAGEMENT

EDGING STRIP can be disposed of with normal household waste after use.

USER INFORMATION

The qualified professional responsible for the preparation, application, or integration of the product into the final application shall provide the user with appropriate information regarding its use, maintenance, and limitations. The user should follow the instructions provided for the correct use of the final application. Contact the qualified professional if adjustments, modifications, or further guidance are required.

STORAGE

- It is advised to store the product in a dark, cool, dry place in the original packaging at the temperature of min. 10°C (50°F) and max. 30°C (86°F).
- Do not lose the labels, They are crucial for product traceability.

GENERAL SAFETY ADVICE

- For external use only. Do not use it on open wounds or in the mouth.
- EDGING STRIP may only be modified or incorporated by qualified professionals.

ADDITIONAL INFORMATION

For additional information, such as distributor contact information, product brochures, safety data sheets, and regulatory information, please visit our website at www.orfit.com

Note:

It is prohibited to make alterations to this text without prior approval from Orfit Industries.

ORFIT INDUSTRIES

Vosveld 9A | B-2110 Wijnegem | Belgium
T (+32) (0)3 326 20 26
welcome@orfit.com

www.orfit.com



Ref. No. 50082

VERSION 4

LAST UPDATE: 13/07/2026

REVISION DATE: 13/07/2028